


DAILY Infusion Lines

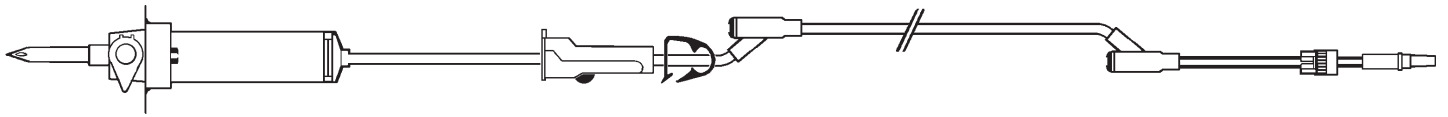


Description

The Daily infusion range are quality one-piece sets offering twin injection ports and set lengths to suit your needs. They are available in a choice of chamber drop sizes.

Daily sets provide precision flow control and feature fully kink-proof end tubing.

Features & Design

- Choice of chamber size (20- or 60-drop/ml)
- Twin Y-port injection sites: easy IV access & low priming
- Precision roller valve and C-clamp for flow control
- Kink-proof end section tubing
- Chambers are filtered and vented
- Fully mobile luer-lock collar allows luer slip attachment
- Convenient port location allows low priming volume i.v. injections as well as i.v. access from outside the kennel
- Clear colour-coding of sets and pouches for identification
- Compatible with standard veterinary infusion pumps

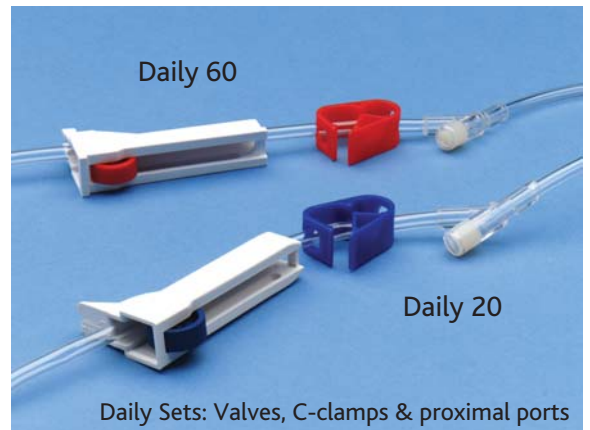

Chamber drop size

Your choice of:

- 20-drops/ml
- 60 drops/ml (micro-dripper)

The 60-drops/ ml chamber makes drop counting easier at low flow rates.

Kink-Proof tubing

End luer connected to a 22G catheter and looped back on itself.

Tube lumen is patent and loop diameter is very small.

Our sets are designed to withstand the rigours of veterinary practice.



Notes: (applying to ALL infusion lines)

1. 'Ports' refers to position (P = Proximal; D = Distal) and type ('Y' = Y-port; 'T' = T-port; '0' = No port) of ports on the set
2. 'Priming' refers to priming volume (ml) from the distal port to the catheter
3. 'Size' is the recommended patient weight range for that set, and is for guidance only. Choose your set based on patient's fluid requirements

Product Summary and Codes

Set	Drops per ml	Length (cm)	Ports P / D	Luers in-line	Priming (ml)	Size	Product Code
Daily 20	20	240	Y / Y	0	0.5	> 10 Kg	I20240YY
Daily 60	60	175	Y / Y	0	0.5	< 20 Kg	I60175YY