


# MODUFLO Infusion Lines

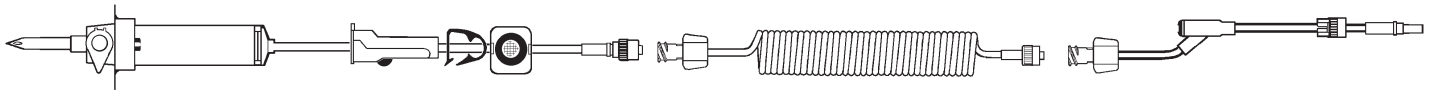

## Description

The ModuFlo range offers the maximum versatility to adapt or alter your infusion set. Twin in-line luer connections allow addition of regulators, valves or extensions, even after the infusion has started. All sets have two injection ports for optimum I/V access and low priming volume, and include a non-kinking end extension.

Tubing on the range has been optimised for accuracy on standard infusion pumps, with the SpiroFlo set being the only spiral line that can be used accurately on infusion pumps. The range can also be used on Baxter pumps.\*

## Features & Design

- Three-piece modular infusion lines
- Twin in-line luer connectors allow you to alter each set to suit every patient, even with the infusion already in progress
- Choice of FOUR sets and three chambers to suit every situation:
  - **20-drop/ml:** Supreme and Spiral (SpiroFlo) sets
  - **60-drop/ml:** Paediatric and Burette sets
- Quality components and colour coding for easy identification
- Kink-proof short end extension as standard on all sets
- Mix & match your components to suit your needs
- Replace tube sections, if required, rather than the whole set
- Add valves, regulators and extensions to improve infusion control
- Compatible with the Optima VS and veterinary infusion pumps

Optima VS

Niki IV

Vet Pro

\* Use on non-dedicated Baxter Flogard models requires inclusion of the Baxter pump extension (EPUMP-BX) at the proximal luer connection



ModuFlo sets are supplied with the rugged short end section as standard



Colour-coded valves, clamps and proximal T port

Burette ModuFlo




150 ml burette body with handle to hang from same hook as fluid bag.

SpiroFlo Standard



Unique spiral tubing gives improved recoil but allows full extensibility. 390 cm long. Fully infusion pump compatible.

Paediatric and Supreme

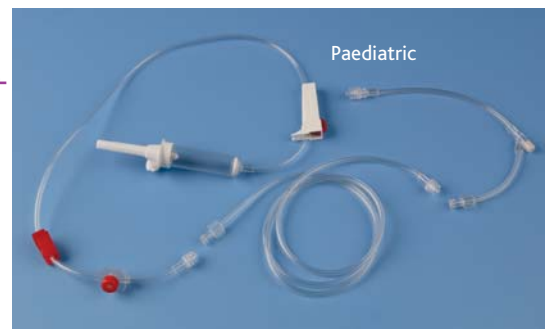

Standard chambers and straight tubing. Both sets offer highly twist-resistant mid-section tubing.

## Product Summary and Codes

Set	Drops per ml	Length (cm)	Ports P / D	Luers in-line	Priming (ml)	Size	Product Code
Paediatric	60	210	T / Y	2	0.5	< 30 Kg	IPM200TY
Supreme	20	280	T / Y	2	0.5	> 10 Kg	ISM260TY
Burette	60	210	T / Y	2	0.5	< 30 Kg	IBM200TY
SpiroFlo	20	390	T / Y	2	0.5	> 10 Kg	ICM20370