

Equine Y-Spiral sets

Design

We offer a choice of TWO Equine Y-Spiral sets to provide maximum versatility with large animal infusions.

Both sets feature twin spike bag access with C-Clamps, fully bonded connections, a quality 10-drop/ml chamber, a T-port injection site and toughened spiral tubing.

Descriptions



Basic Y-Spiral set (IEQYB-01)

- One-piece 10-drop/ml Y-spiral infusion line
- Specially toughened PVC spiral tube section
- Responsive Roller valve and C-clamp for guaranteed flow control
- Bore diameter: 4.8mm
- Tie wraps can be supplied for securing spikes and line if required

Fulmac Modular Y-Spiral set (IEQYPUR2)

- One-piece 10-drop/ml Y-spiral infusion line
- Specially toughened Polyurethane spiral tube section for recoil
- Responsive Roller valve and C-clamp for guaranteed flow control
- Twin in-line luer connections for modular functionality (see ModuFlo sets)
 - ⇒ Add valves, connectors or extension lines
 - ⇒ Add swivel connector to prevent line tangling
 - ⇒ Simple disconnection from patient
 - ⇒ Convert to Wide Bore Y-set using the Fulmac extension (ESWB2200)
 - ⇒ Set includes it's own end extension with C-Clamp attached
- Ancillary line included for secondary infusions or drug administration
- Bore diameter: 4.3mm
- Supplied with tie wraps for securing spikes and line if required

* Ask about our transfer sets for prolonged fluid infusions on these twin-spike sets *

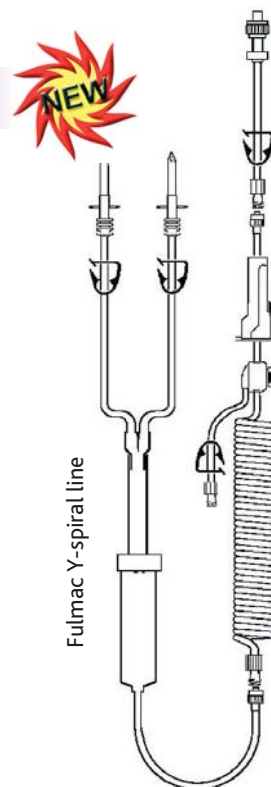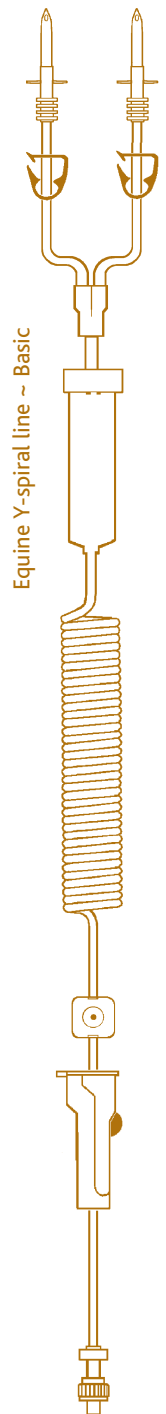

Product Summary and Codes

Set	Length (cm)	Ports P / D	Luers (in-line)	Flow Rate * (L/hr)	Product Code
Basic Equine Y-spiral	550	0 / T	0	25-30	IEQYB-01
Fulmac Modular Y-spiral set	570	0 / T	2	20-25	IEQYPUR2

* Typical flow rate at 100 cm drop with no catheter connected