

END Extensions

Design

End extensions allow quick and easy disconnection from the set without touching the catheter connection. They can all be secured above the catheter site, reducing the risk of catheter movement and making patient mobility so much easier.

All extensions are designed to give very low injection priming volumes but be robust enough for even the most vigorous patient. They allow the addition of valves and connectors, even whilst infusions are running, to make managing the infusion easier and safer.


We offer a choice of FIVE designs and cover three injection port types: Standard 'Y', needle-free 'Y' and 'T'.

Description - general


- Low or ultra-low priming volumes on all extensions
- Choice of FIVE extensions and THREE port types
- Allows placement of valves & connectors to protect your patient and make the infusion easier
- Improves patient mobility, saves you time and reduces risk of displacing catheter when handling patient
- Ideal for easy IV access even when not running fluids (eg frequent IV drugs, anaesthesia or blood sampling)
- ALL ModuFlo sets use EEKP4125 as a standard end

**** NOTE: When disconnecting set from end extension, the female luer must be capped to maintain sterility and prevent air ingress. Using our Needle-Free valve (FNV-MF-H or FNV-MULTI) is the easiest and safest way of doing this. Please refer to the 'Making the most of ModuFlo' guide for further details. ****

Mini-bore End Extensions

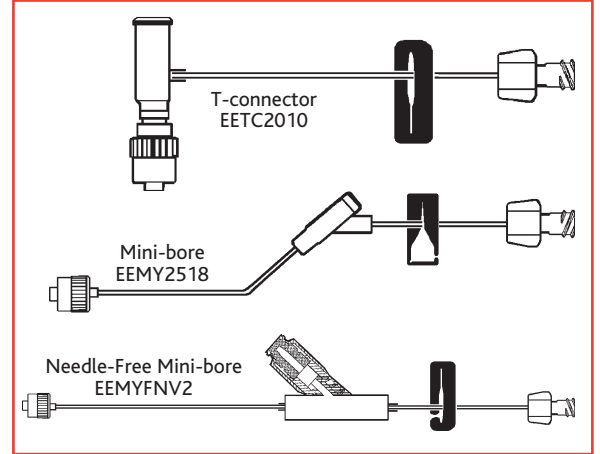

Kink-proof End Extensions

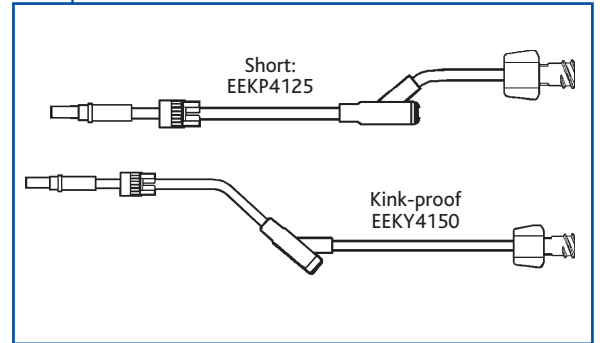

NON-Kinking Extensions

Features & Descriptions


- Male luer has fully mobile collar for push-fit or luer-lock
- Extremely tough tubing but good pliancy, allowing use in smaller patients too
- Clinically kink-proof between injection port and male luer to ensure the fluid path remains open at all times
- Standard non-latex 'Y' port injection site 20-25 cm from end

★ SHORT End extension (EEKP4125)

- ⇒ Standard end extension on all ModuFlo sets
- ⇒ Kink-proof tubing above and below injection port
- ⇒ Compact 25 cm length

★ Kink-Proof End extension (EEKY4150)

- ⇒ Ideal for large animal practice when high flow rates are not critical
- ⇒ Kink-proof tubing below injection port, kink-resistant above
- ⇒ Longer 50 cm extension

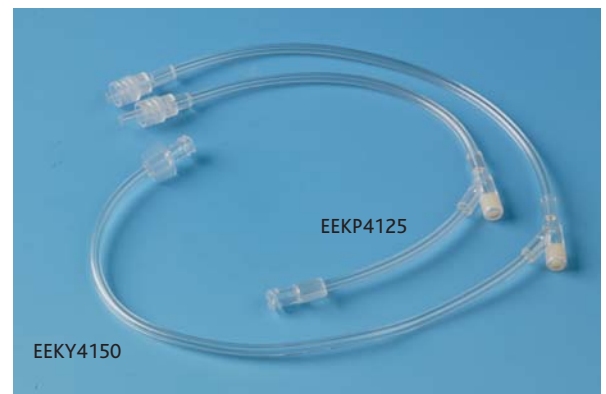

Product Summary and Codes

Extension	Length (cm)	Priming (ml)	Ports (type)	ModuFlo end-section	Product Code
Short End	25	0.4	Y	ALL	EEKP4125
Kink-proof End	50	0.45	Y	-	EEKY4150