

## Standard Drain Connectors

### Features

- Connectors to fit most standard drain tubing
- Your choice of MALE or FEMALE luer connection
- Supplied as single sterile units

### Details

- Fits tube sizes 12-28 Fr
- Female connector allows attachment to standard extensions
- Male connector allows reversed attachment to suction pumps etc



## Needle-Free Drain Connections

### Features

- Fully assembled ready to fit to drain
- Choice of three configurations to suit most situations and sizes
- Fitted with a needle-free valve for simple access and closure
- NO risk of ever leaving a tap handle open!
- Valve is swabable, so does not need capping between uses
- C-clamp fitted on connector tubing for closure if changing valve
- Additional C-clamp included for use on chest drain
- Can add valves (eg Heimlich) etc as needle-free valve is removable
- Your complete drain closure set in a single connector!

This drain connector is closed with the needle-free valve (FNV-MF-H).

It gives the most secure closure whilst allowing the easiest access for drainage.

Connect male luer to the valve to start drainage/aspiration.



### Needle-Free Drain Connector (DR-FNV-1)


- Fits same tube sizes as standard drain connectors above
- Compact & light weight
- Suits low to moderate flow rates; pneumothorax and transudate
- Best rescue/replacement connector if any connector part is damaged

### Wide Bore Needle-Free Connector (DR-FNV-2)

- High flow rates, suits thick exudates
- Use on 16 - 24 Fr drain sizes only

### Needle-Free Drain Adaptor (DR-FNV-3)

- Universal connector for any drain with standard connector
- Ensures you always have a connector to fit your drain

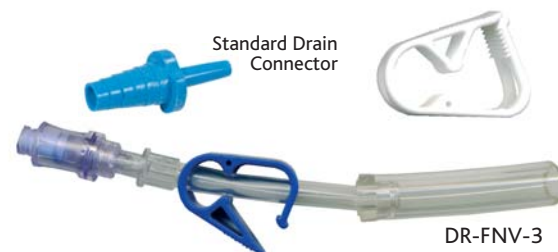

## Product Summary and Codes

| Drain Connector                 | Luer   | Drain Size (Fr)      | Needle Free? | C-Clamp Fits (Fr): | Product Code |
|---------------------------------|--------|----------------------|--------------|--------------------|--------------|
| Standard; Female luer           | Female | 12 - 28              | No           | n/a                | C-DR-FLL     |
| Standard; Male luer             | Male   | 12 - 28              | No           | n/a                | C-DR-MLL     |
| Needle-free drain connector     | Female | 12 - 28              | YES          | 12-18              | DR-FNV-1     |
| Wide Bore Needle-free connector | Female | 16 - 24              | YES          | 16 - 28            | DR-FNV-2     |
| Needle-free drain adaptor       | Female | Any (with connector) | YES          | 16 - 28            | DR-FNV-3     |