


## Oxygen Connector Extension

- Luer-connected extension to connect oxygen supply to nasal tube/prongs
- Generous 210 cm length; can be altered to suit patient by changing extension
- Light weight and fully adaptable whilst in use
- Excellent flow rates
- **FITS** to oxygen supply as follows:
  - ⇒ **Main adaptor tubing:** - Direct to oxygen spigot
  - 8 - 8.5 mm endotracheal tube connector
  - ⇒ **Secondary tubing:** - Cut to fit 5 mm endotracheal tube connector
  - ⇒ **Extension tubing:** - Cut to fit 3.5 mm endotracheal tube connector
- **CONNECTS** to delivery device as follows:
  - ⇒ **Cone/Funnel devices:** - Use Cone Connector (C-CONE-1) - see below
  - ⇒ **Luer-connected devices:** - Use male luer on extension tubing
  - ⇒ **10 Fr nasal catheter:** - Use male luer of extension tubing (on catheter NO201110)

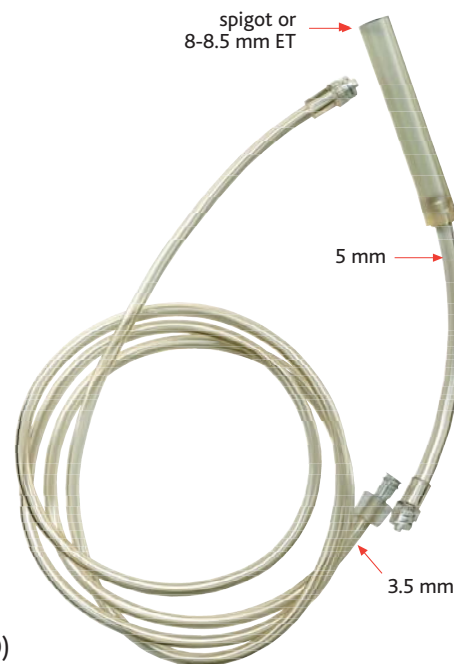

## Product Summary & Code

Product	Length (cm)	Pack Size	Product Code
Oxygen Connector Extension	210	1	NOE2100L

## Cone Connector/Adaptor


- Connects Oxygen Extension (NOE2100L) onto most cone or funnel delivery sets
- We recommend bonding the adaptor to the delivery device for security. The oxygen line can still be altered via the standard luer connections
- Fits directly to:
  - ⇒ **Nasal Oxygen Prong sets** (Adult & Paediatric) : - NP-2008A
  - NP-2008P
  - ⇒ **Nasal Oxygen Catheter sets** (6 & 8 Fr) : - NO201106
  - NO201108
- Can also be used in reverse direction to fit drain extensions up to 24 Fr bore size



Cone connector (C-CONE-1)



Cone connector with nasal oxygen prong funnel end

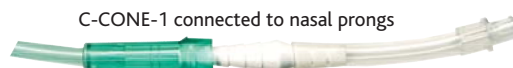

C-CONE-1 connected to nasal prongs

## Product Summary & Code

Product	Length (cm)	Pack Size	Product Code
Cone Connector	9	1	C-CONE-1